

NURIKON Aligner Pliers are specially designed to prevent distortion of the aligner without the need to heat the plier. These pliers allow greater torque, rotational control and increased retention.

Family Name: Orthodontic Pliers

Intended Use: The plier is used for applying a hole or shallow indentation on the orthodontic Aligner. Its use is limited to dental orthodontic treatments.

Thermal and Non Thermal Pliers.

Non Thermal Pliers

Non Thermal Pliers

The Hole Punch Plier

NK-261 Punch Size 6mm

NK-261 A Punch Size 5mm

NK-261 B Punch Size 4mm

The Hole Punch Plier is used to cut half moon shaped hole at the gingival margin of the clear aligners. These half circles permit the addition of bonded buttons, bonded orthodontic tubes or brackets with associated hooks or as simply used to relieve impingement of plastic on soft tissue.



THE HORIZONTAL / Torque Plier

NK-262

The Horizontal Aligner Plier creates a horizontal indentation, for individual root torque. It also enhances retention of either clear aligners or retainers by placing an indent into a tooth's undercut or below a bonded attachment.



THE VERTICAL / Rotation Plier

NK-263

The Vertical Aligner Plier creates a vertical indentation for individual tooth rotation. The Vertical can also be used to produce an indent to produce a mild "nudge" to enhance resolution of in-and-out discrepancies.



PEARL DOT Smooth 1mm / Retention Plier

NK-264

The Pearl Dot Smooth Aligner Plier creates a circular indentation, 1mm in diameter, for increased retention in retainer and aligner system.



PEARL DOT Smooth 2mm / Retention Plier NK-264-A

TEAR DROP/ Key Hole Cutter

NK-265

Tear Drop plier cuts tear-shaped notch or hook at gingival margin of clear aligner. It provides a reservoir to hold the elastics in place as the aligner is seated. This permits easier manipulation of elastics by the patient.



Dimple Erasing Plier

NK-277

This Plier is used to remove protrusions and dimples on the aligners. Simply place the beak on the dimple and squeeze closed to flatten.



Thermal Pliers



**Circular Thermal
Forming Plier
NK-272 C**



NK-274 E



**Multi Indent
Forming Plier
NK-271 D**

Bubble Forming



NK-276 B

Under Cut Enhancer



NK-270 A





PUNCH PLIER
NK-266 A 2mm
NK-266 B 3mm
NK-266 C 4mm



Slot PUNCH Plier
NK-269




Serrated Trimming Scissors - NK-279

- Serrated design for controlled trimming without slippage
- Curved tip cuts easily and precisely, even in hard to reach lingual areas
- 3 1/2" curved serrated trimming scissors

NK-278

These lightweight scissors are used for removing excess thermoformed material, and for the smooth trimming and shaping of aligners. Their streamline design, graduated curve and tapered beak makes for easy use - even in tight areas.

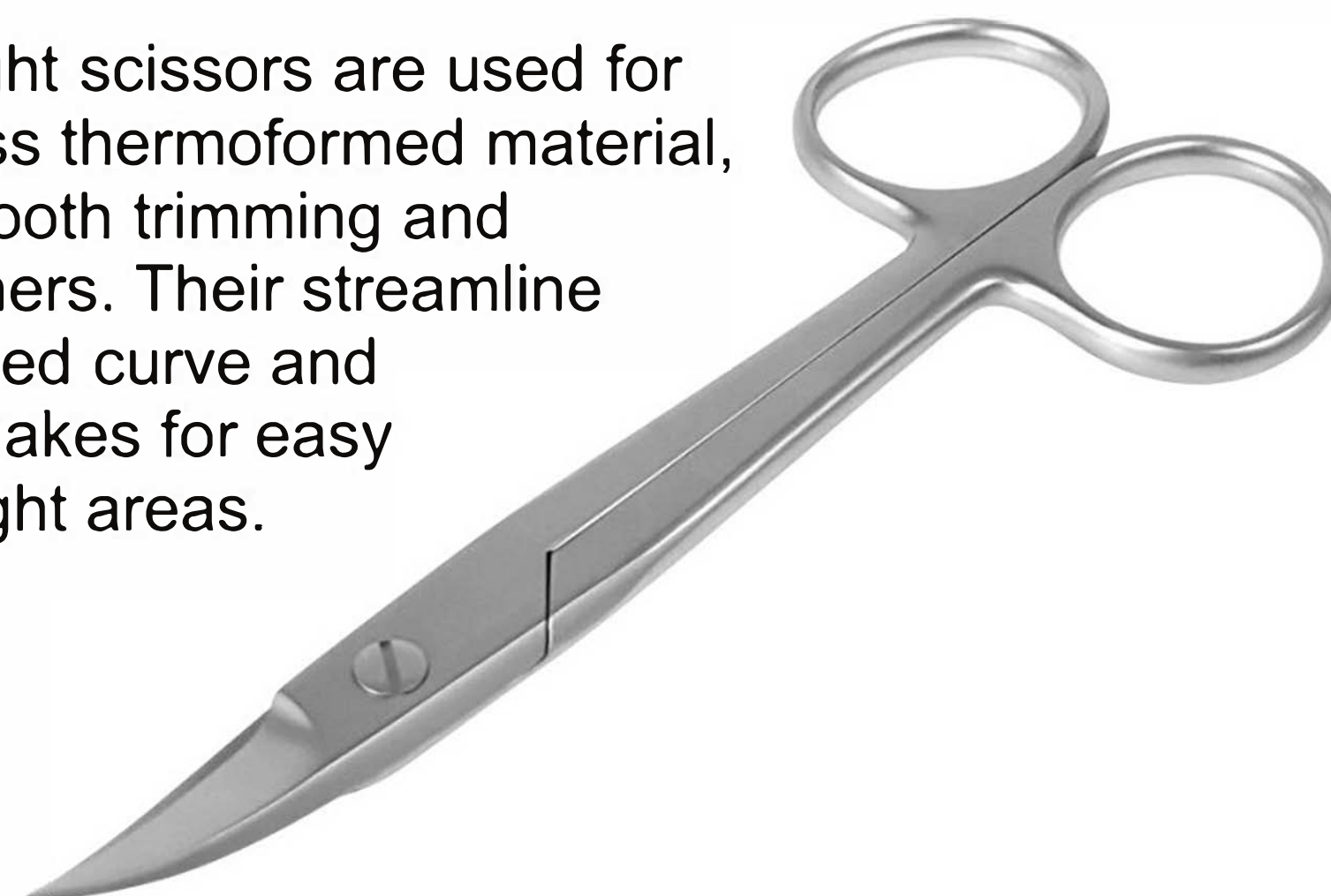

Heavy Aligner Scissors, Serrated

NK-404





NK-333-A
Height Gauge Adjustable



NK-333-C
IPR
(Interproximal
reduction)

Gauge Set Used during Interproximal Reduction (IPR)
to measure and confirm space.

Stainless steel, 1 per pack Six Gauges Included:

10 mm,



20 mm,



25 mm,



30 mm,



40 mm,



50 mm



NK-401
Instrument
for opening
and losing
the brackets



NK-403
Tweezer



NK-402
Self Ligating Tool

This is Cfast's self ligating too to be used with the self ligating brackets.

**RAPID PALATAL
EXPANDER
(RPE) Removing Plier
(Acrylic Debonding Plier)**

NK-533

Designed to quickly and
easily remove bonded
orthodontic acrylic
appliances,
such as RPE's and
Herbst.
The plier has a hinged
platform on one side and
a debonding blade on the other.

